

REACH B-LI RUGGED

OUR TOUGHEST HEARING AID NOW WITH MULTI-VOICE FOCUS



- › Enabling users to be fully involved in group conversations with Multi-Voice Focus
- › With My Voice 2.0 the users own voice stays comfortable with a dedicated processing separated from the target speech sounds
- › Handsfree phone calls via Headset Mode with Call Control when connected via MFi or LE Audio
- › Rechargeable solution for a full days use:
 - › Rugged delivers up to 32h runtime (incl. 5h streaming)
 - › Rugged delivers up to 36h runtime without streaming
- › Lifeproof to withstand everyday challenges
- › AI with Rexton Assist helps you to seamlessly customize your hearing aids


Rugged offers a uniquely tough build:

- › Water resistant construction
- › Sweat proof
- › Drop proof
- › Scratch resistant
- › IP68 certified
- › ATEX certified*

Colors

Dark Champagne Granite	Granite Black	Grey Granite	Rose Gold Rose Gold	Sandy Brown Granite	Silver Silver	Beige Beige	Black Black
--------------------------	-----------------	----------------	-----------------------	-----------------------	-----------------	---------------	---------------

Fitting range is based on 2 ccm coupler:

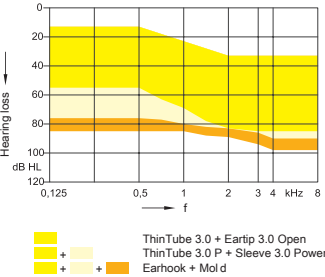

- Earhook
- 60 dB / 133 dB SPL (2 ccm coupler)
 - 67 dB / 138 dB SPL (Ear simulator)

TOUGH DESIGN WITH A SUPER-COMPACT SIZE

Our toughest hearing aid is now available within our platform Reach. Reach Rugged will deliver excellent sound processing in a small and secure package. It's structure has an airtight receiver, sealed charging contacts, reinforced receiver suspension, thickened housing, and sealed parting lines.



Features / Performance Level		80	60	40	
Channels/Controls/Programs		48/20/6	36/18/6	32/16/6	
LED		-	-	-	
T-coil		•	•	•	
Bluetooth		•	•	•	
Tinnitus Sound function		•	•	•	
Tinnitus Notch function		•	•	•	
Wireless sync		•	•	•	
IP68		•	•	•	
XPhone		•	•	•	
Technology Highlights					
Multi-Voice Focus		•	•	•	
Focus Beams		4	2	1	
Voice Stabilizer		•	•	•	
Speech Preservation Technology		•	•	•	
My Voice 2.0		•	•	•	
Motion sensor for scene analyzis		•	•	•	
Directionality					
Automatic adaptive		•	•	•	
iOmni Sound Locator		•	•	•	
Auto iFocus 360		•	-	-	
iFocus 360		•	•	-	
Auto Focus 360		-	•	-	
Auto Stereo iLock		•	•	•	
OneMic Directionality		-	-	-	
Noise Reduction					
Noise Management		•	•	•	
Sound Smoothing		•	•	•	
Wind Noise Reduction		•	•	•	
Auto Echo Reducer		•	-	-	
Reverb Reducer		•	•	-	
Sound Quality & Control					
Bandwidth extension		•	-	-	
Bandwidth compression		•	•	•	
Feedback Preventer		•	•	•	
Music Enhancer (presets)		3	3	1	
Automatic Acclimatization		•	•	•	
Audio Streaming					
Bluetooth Classic		-	-	-	
LE Audio		•	•	•	
MFi		•	•	•	
ASHA		•	•	•	
Headset Mode					
Bluetooth Classic		-	-	-	
LE Audio		•	•	•	
MFi		•	•	•	
TeleCare					
Full live remote tuning		•	•	•	
Basic portal remote tuning		•	•	•	
Accessories & Apps		Chargers		CROS	
Smart Key	•	Charging Station	•	Reach CROS R-Li T BC	-
TV Sound	•	Travel Charger	-	Reach CROS R-Li	•
Smart Transmitter 2,4	•	Charging+ Station	-	Reach CROS SR	-
Smart Mic	•	IF CIC Travel Charger	-	Reach CROS iX-CIC Li	-
Rexton App	•	Universal Cable Charger	-		
Programming					
Noahlink Wireless	•				
ConnexxAir	-				
HI-PRO	-				